

# PINNACLE SL300

Bringing You Peace of Mind



## DIFFERENT BY DESIGN

The Pinnacle Stair Lift has an established reputation as an industry leader for its innovative design, solid performance and reliability. It's patented drive system serves as the cornerstone that makes this stair lift different from all others on the market.

- Lubrication free, keeps you, your stairs and even your pets free from a mess
- Low friction system provides reduced cost of ownership: no maintenance required and lower energy bills.
- Runs on batteries with efficient motor resulting in extended performance during power loss
- The over-sized worm gear offers greater stability
- Symmetrical design provides industry's most compact folded width



## CLEAR WALKWAYS

Harmar's exclusive automatic and mechanical folding rail removes obstructions and trip hazards at the lower landing area. Its operation is automatic, the user does not need to stop or have someone else assist. And since it's mechanical, there are no additional switches or controls.

# DESIGNED WITH YOU IN MIND

## Pinnacle SL300 Stair Lift




Ergonomic control switch helps those with arthritis or limited hand dexterity



Oversized footrest provides more comfort and security



Nylon-polymer gear with helical worm drive features 6 points of contact for ultimate stability

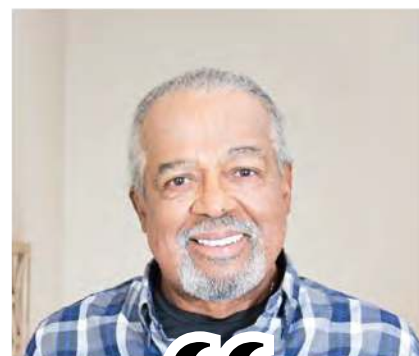

*My Pinnacle gives me the confidence to go up and down my stairs more safely. It operates so smoothly and is very quiet. I now feel more comfortable in my own home and less dependent upon others. My children have peace of mind knowing I am safer with my Pinnacle.*

# 40+

Cycles During Power Loss

# 300

Pound Capacity

# ZERO

Maintenance Required

# 11"

Folding Width from Wall

### SL300 PRODUCT SPECIFICATIONS

Capacity	300 lbs.
Speed	20 fpm
Folded Width	11"
Drive System	Helical Worm Gear with Nylon Polymer Rack
Power	24V DC (Battery), Charger plugs into standard AC outlet
Rail Length	16' - 0" standard
Safety Feature	Obstruction sensors, 90° seat swivel interlock, Seatbelt
Warranty	2 Year Parts, 10 Year Worm Gear / Rack, 1 Year Batteries
Option	Retractable seat belt, Key lock, Extra rail

### SL300FR PRODUCT SPECIFICATIONS

Folding Rail	Automatic, mechanical & continuous motion
Folded Width	13"

### COMPACT

The Pinnacle SL300 chassis folds up neatly, occupying a mere 11" of width from your staircase wall, allowing others to pass by. The rail mounts only 1" from the wall and is only 3" tall helping to stay out of sight.

